



DISPOSABLE NOVA CLIP APPLIER

Nova Surgery

ENDO NOVA CLIP APPLIER



10 MM SHAFT



ML CLIP SIZE



ROTATION KNOB



TRIGGER



COUNTER



HANDLE




DESCRIPTION


The NOVA Clip Applier is a sterile, single-use endoscopic tool designed for reliable ligation in surgical procedures. Its enhanced jaws and refined design minimize clip scissoring and splitting, offering improved grip and secure clip placement on vessels. Pre-loaded with 20 titanium clips that advance automatically with each use, it fits standard trocars. The shaft is made from a low-glare material to reduce reflections, ensuring better visibility during surgery.



IMPROVED JAW SECURITY
Features jaws stronger than other brands, ensuring unmatched clip formation

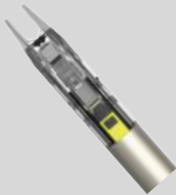


REDUCED VESSEL DAMAGE
visual reference of the quantity of clips remaining with last clip lockout feature

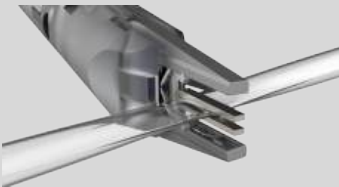

ENHANCED CLIP STABILITY
Reduced clip scissoring by ensuring superior vessel clip formation compared to the previous model.

PRODUCT INSIGHTS

- Last clip lockout with indicator bar
- Anti-clip drop ratchet mechanism
- 360° rotation knob/shaft — designed for easy rotation and one-handed use
- The rigid, low-glare shaft minimizes reflection and movement for precision
- Ergonomic handle for easy control
- Inner-crossed pattern of the clip offers a firm hold on the vessel

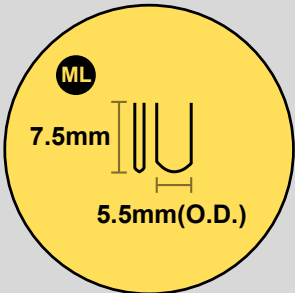

ORDER DATA

Item Code	Description	Diameter	Working Length	No of Tack	Unit / Box
YBJK- L	NOVA Clip Applier	10 mm	360 mm	20	1 pc



Nova Energy

Nova Surgery

 No.51 Shuishan Road, Zhonglou Economic Development District,
Changzhou, Jiangsu, China. Zip code: 213023

 +86 1816 8812 599

 Info@xnystaplers.com

 www.xnymedical.com

